



LOVE MY SKIN



REJUVENATION BEAUTY WAND

LMSRBW

Introduction


Thank you for choosing the Love My Skin Rejuvenation Beauty Wand to achieve more radiant and flawless skin.

Please read the instructions before you turn on this device.

This device uses LED light therapy and EMS Electrical Muscle Stimulation to gently activate the skin cells, stimulate rejuvenation of collagen, improve skin elasticity, soften wrinkles and leave the skin looking smooth and healthy.

This painless treatment can plump the skin and treats mild to moderate acne.

The LED light works by sending light into the deeper layers of the skin and can minimize fine lines, treat sun-damage and stretch marks and reduce redness after more aggressive IPL or laser treatments.

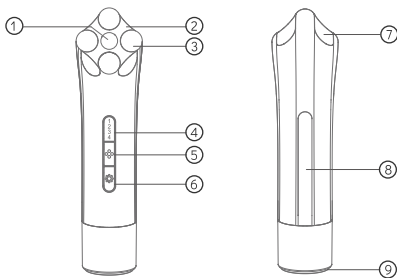

① Massage Head

④ EMS Level Indicator

⑦ Charging Indicator

② LED Lighting

⑤ EMS ON/OFF

⑧ Conduction Panel

③ Massage Tips

⑥ Phototherapy ON/OFF

⑨ USB Charging Port

WARNING: NO MODIFICATION OF THIS EQUIPMENT IS ALLOWED.

Settings

This device has three modes: EMS Mode, Phototherapy Mode, and EMS+ Phototherapy Mode.

EMS Mode = anti-aging electrical muscle stimulation to tone.

Phototherapy mode = anti-aging light therapy, to reduce acne and promote clearer skin.

EMS + Phototherapy = to promote younger, toned, clearer, smoother, healthier, flawless skin that glows.

- Simply press button 5 once to activate the EMS Mode, and use the same button to select different intensities. Press and hold the button 5 to deactivate the EMS Mode. Note: Conductive Gel or non-oil based serum should be applied before using the EMS Mode. It is recommended to use level 1 for the first time and then enhance the level as needed.
- Simply press button 6 once to activate the Phototherapy Mode, and use the same button to select different light therapy. Press and hold the button 6 to deactivate.

Please refer the table below to choose the right light therapy for you.

Light Therapy	Purpose
Blue Light	Soothes and diminishes skin inflammation
Yellow Light	Treats hyperpigmentation
Red Light	Promotes collagen regeneration
Combination	Combination light with multiple benefits

Note: Avoid direct eye contact with the lighting source while using the Phototherapy Mode.

- For maximum result, turn on both EMS + Phototherapy Modes, and choose the right intensity and lighting therapy according to your need.


How To Use

1. Cleanse your face to remove any makeup, oil and dirt.
2. Please change step 2 to "2. Apply your favourite oil free skin care serum and or cream evenly on your face, and ensure you keep your face moisturised whilst using our rejuvenation beauty wand.

Warning: If you decide not to use conductive gel or serum whilst in EMS mode you may feel a slight tingle in your skin, we strongly recommend you keep your face moisturised with your favourite skin care serums and or creams whilst using this device to enjoy this treatment and improve the EMS technology benefits.

3. Turn on the device and choose the mode. Starting from the forehead, gently press the massage head against the treated area.


Move the device in circular motion as refer to the diagram below:

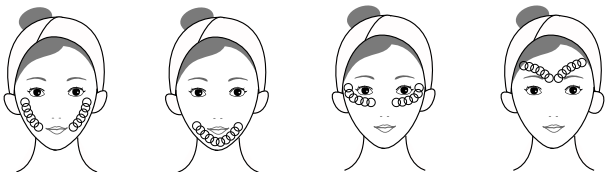

4. Turn off the device and clean the massage head after each use.

Note: It is recommend to use the device 2-3 time a week, and up to 10 minutes each treatment.

TIP Make sure the massage head is fully in contact with the treated area. Please refer to the diagram on the right:

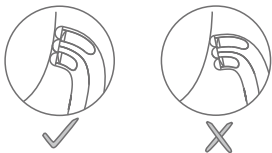

Charging Instructions

Please note charging can take up to 3 hours.
The full charge should last for over 2 months of use.

1. Please charge the device before first use, charging time is approximately 3 hours.
2. When charging, the red light will be on. When the device is fully charged, the light indicator will turn green.
3. The device will beep three times and then deactivate automatically when it runs low on battery.

⚠ WARNING: Before charging, please make sure the charging cable and device are away from water and extreme heat.

Warranty

1-Year Warranty

One-year warranty is provided on this device from the original date of purchase, warranty is valid when product is used according to the instructions.

The warranty covers working parts that affect the function and does not cover cosmetic deterioration caused by wear and tear or any damage caused by accident, misuse or neglect.

Any attempt to disassemble this device will void the warranty.
If you discover a defect and notify us during the warranty period, we will, at our discretion, replace the device.

Claims under warranty must return the faulty device with proof of receipt within the warranty period.

Warnings

1. Do not use this device near water or unsafe surroundings.
2. Do not use this device if you have sunburn, broken or irritated skin, swelling or have had skin related surgery.
3. Do not use this device if the device is damaged.
4. Do not disassemble or modify this device.
5. Avoid any contact with your eyes open.
6. Keep out of reach of children.
7. Discontinue use if you experience any pain or discomfort.
8. Do not use this device on any infected areas of skin such as rashes or open wounds.
9. Do not use this device during pregnancy.
10. Do not use this product near flammable substances such as alcohol, oils and aerosol sprays etc.
11. Do not use if you are sensitive to light.

Device Specification

Model	LMSRBW
Rated operational voltage	5V
Rated operating power	<5W
Rated operational current	<1A
Power supply	USB charge
Weight	115g
What's in the box	Rejuvenation Beauty Wand Device, USB Cable, User Manual

Please feel free to contact us at:



+44 (0) 1483 937967



hello@love-my-skin.co.uk



Love My Skin (Kinetik Medical Devices)
Unit 3, Perrywood Business Park
Honeycrook Lane, Salfords Redhill
RH1 5DZ. UK



LOVE MY SKIN



Symbol for "Manufacture Date".



Symbol for "Environment Protection" – Waste electrical products should not be disposed of with household waste. Please recycle where facilities exist. Check with your local Authority or retailer for recycling advice.



Symbol for "Complies with EU requirements".



Symbol for "Complies with UK requirements".



Follow this instruction manual for your safety.



Model Reference.