

SKYN CARE.co

BLACKHEAD REMOVER

INSTRUCTION MANUAL



LMSBHR1
LMSBHR2

PRODUCT INTRODUCTION

Microcrystalline Firming Blackhead Remover is a multifunctional beauty instrument which integrates the functions of microdermabrasion, firming skin, cleansing pores, blackheads and acne. The vacuum aspirator utilizes more than 100,000 natural mineral micro-drilling particles to remove the aging skin and the dirt in the pores reduce facial fine lines, which is a non-invasive and non-irritating technique, uses the diamond thickness on the suction bar to control microdermabrasion strength. At the same time, 4 different shape probes have 4 different functions, which combines various functions such as microdermabrasion, firming and clean pores.

WORKING PRINCIPLE

Vacuum vapour-pressure type slimming face technology

1. With a vacuum suction system, it gently massages the skin, promotes blood and lymph circulation, promotes metabolism, and makes the dermal tissue get enough nutrients to make the skin more firm and shiny;
2. Enhances the absorption of your favorite beauty products so they can penetrate the skin;
3. Promotes collagen reproduction, enhancing the skin to look young and firm;
4. Promotes skin cell regeneration;
5. Improves skin microcirculation and promotes melanin metabolism, thereby diminishing melanin spots on the skin.

Microdermabrasion probe

The natural mineral microcrystalline diamond particles on the microdermabrasion probe gently remove the stratum corneum to make the skin smoother and more radiant. It gently pushes out the rough surface and debris, and has the function of suction. It can suck out the dirt on the skin, promote the blood circulation while removing the skin old dead cell, help skin cell renewal, keeps the skin smooth and rejuvenates.

PRODUCT DETAILS

Name: Blackhead Remover
Power input: 5V•1A
Power: SW
Size: 140.35mm
Body weight: 67.6g

PRODUCT FUNCTION

Four replaceable heads with different functions:



Microcrystalline skin head

A microdermabrasion probe with microcrystalline diamond particles. The surface area of the crystallites is large, and the surface adsorption performance and activity are quite prominent. In the case of not harming the skin, it is better to massage the skin to achieve the exfoliating effect, so that the skin is rejuvenated!



Blackhead removal & skin firming head

The large round mouth tip has the characteristics of strong suction. It has a remarkable effect on fat blackheads and so on. The unique vacuum pressure suction mechanism is used for micro-lifting and massage on the skin, deeply affecting the bottom layer of the skin, stimulate microcirculation operation, promote blood and lymph circulation, and recreating firm skin!



Clean pores head

The small round hole probe has a weak suction force. Aim the round holes at blackheads and acne to remove dirt, blackheads and acne. If the acne is severely broken and redness, it is not recommended. If the blackhead is stubborn, it is recommended to use a large round hole probe. Suitable for thin and sensitive skin.



Reduce fine lines head

The special elliptical hole design is designed for micro lifting massage of unique parts of the face such as the corners of the eyes, the nose and the corners of the mouth, which are difficult to care for, and deeply acts on the bottom layer of the skin to stimulate microcirculation and prevent deep wrinkles.

A SAFETY WARNING

1. The firming function has the function of facial V shape lifting, and skin with injuries and scars is not recommended.
2. Dry skin is recommended to be used once a week. Oily mixed skin can be used twice a week for less than 5 minutes each time, but it is also depending on the skin condition.
3. Please do not press hard on the skin, the instrument itself has a strong suction.
4. After the end of the nursing, the skin will have some irritation or redness due to the accelerated microcirculation of skin. This is a normal phenomenon. After stopping use, it will automatically disappear after about 30 minutes.
5. Make sure you are moisturizing and using sunscreen when you have finished using the device.
6. Do not use scrubs or other exfoliating products before and after treatment.

HOW TO USE

FACIAL MICRODERMABRASION AND FIRING SKIN FUNCTION:

1. Make sure you have a clean and dry face. Choose which sized beauty head you want to use and attach to the device.
2. Starting from the jaw, move the beauty instrument from the bottom up. Place the device on the cheek and move from the inside out.
3. Place the device on the forehead and move from the middle to the sides then place the Blackhead Remover in the T zone and move from the bottom to the top of the face.

HOW TO REMOVE BLACKHEADS

1. Replace with a large round hole or a small round hole beauty head.
2. Use the device on the blackheads and acne. Move the device around your face as needed. .

Technique: When you suck the blackheads move the device around slowly around the targeted area but please don't move back and fourth.


REDUCE FINE LINES

1. Insert the oval hole beauty head.
2. Place it in the corners of the eyes, nose, mouth, etc. for care, adjust the machine position as needed.



Technique: You can use it to suck the skin parts that need to be treated, and then quickly pull it out, and repeat this operation many times. This operation can act on the bottom layer of the skin to stimulate microcirculation, promote blood circulation and increase skin elasticity and prevent wrinkles.

PRODUCT DETAIL

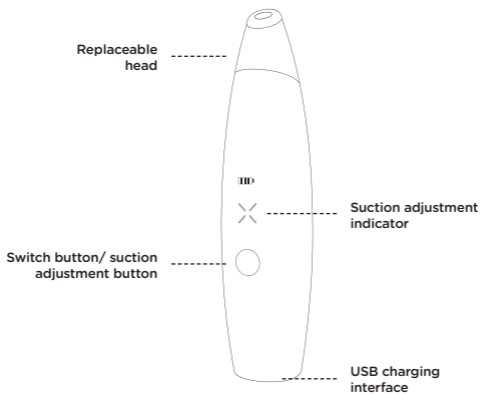

MODE DISPLAY

Display status:
Mode 1



Display status:
Mode 2



Display status:
Mode 3



Display status:
Mode 4



POWER DISPLAY

When the battery is low, all the symbols flash together.

0.2 SECOND INTERVAL



When charging, the symbol is displayed in a clockwise rotation cycle.
When the charging is completed, the symbol is always on.

WHERE TO USE THE 4 SUCTION HEADS

MICROCRYSTALLINE HEAD:

Exfoliating effect for rejuvenating skin

BLACKHEAD REMOVAL HEAD:

Large round head for strong suction

PORE PERFECTION HEAD:

Delicate blackhead removal for sensitive skin

SPECIAL POINTED HEAD:

Designed for Microliftin massage suction for the corners of the nose and mouth.

PRODUCT CONFIGURATION



①



②



③



④



⑤

1. Device x1
2. Suction Head x4
3. Sponge x5 + Silicone ring x2

4. Charging cable x1
5. Instruction manual x1

PRECAUTIONS

1. Ensure your hands are dry before use and do not let water come into contact with the device as it could damage it.
2. Make sure the filter sponge is properly inserted into the filter cup before handling and make sure the sponge is dry. Make sure water or other liquid will not be sucked into the air pump of the body, causing blockage, the device will not work properly.
3. Do not wash the product with water.
4. This product is not intended for use by persons with disability or lack of experience and knowledge, (including children). There must be a guardian who supervises and guides the use of this product and is responsible for how the device is used.
5. This product is not suitable for cosmetic surgery, rhinoplasty and other parts of the skin cosmetic surgery, including pregnant women and heart patients.



Model Reference



Manufacture Date



Complies with EU requirements



“Environment Protection” – Waste electrical products should not be disposed of with household waste. Please recycle where facilities exist. Check with your local Authority or retailer for recycling advice.

