




LOVE MY SKIN

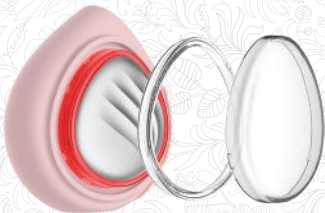

MIRACLE FACE MASK DEVICE

LMSSM

Introduction

Thank you for choosing the Love My Skin Miracle Face Mask Device.


This is your ultimate at home pamper session for an accelerated mask treatment in the comfort of your own home, providing optimal skin care results in minimal time.

The vibrating Miracle Face Mask Device is made from hygienic silicone and innovative photon thermo-therapy technology, gradually warming up to disperse and deeply penetrate the goodness of your favourite facial sheet mask ingredients into the skin, combined with yellow, red and blue LED light which helps to target a variety of skin care concerns, promoting a glowing complexion.

Easy to use, you apply over your choice of facial sheet mask or secure a mask into the device with the detachable ring to commence this next generation treatment of facial indulgence.

Product Diagram

- ① Transparent protective cover
- ② Silicone Body
- ③ Massaging Plate
- ④ Transparent fixed ring
- ⑤ Ring grooves
- ⑥ ON/OFF/Mode Button
- ⑦ DC Port
- ⑧ DC charging cable
- ⑨ Pouch

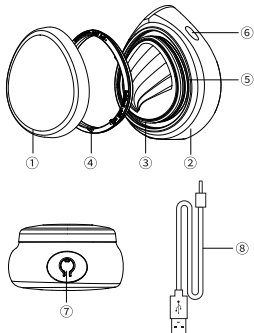

Note: This product comes with built-in lithium and ONLY be charged with USB cable(DC port)

Settings

This Miracle Face Mask Device uses yellow, red and blue LED light.

Yellow LED promotes cell growth, and effectively repairs the underlying skin, we recommend you apply with a hydrating mask.

Red LED light helps to stimulate collagen regeneration and brighten the skin tone, we recommend you apply with a rejuvenation mask.

Blue LED sterilizes and mitigates inflamed acne, stabilizes overactive oil production, we recommend you apply with a soothing mask.

The thermal care helps the skin to open the pores whilst the vibration accelerates absorption of the mask ingredients, effectively shortening the traditional facial mask skin care time.

No.	Item	Indications
1	ON	Quick press the ON/OFF button where indicated on the device diagram. After 1 beep, 1/3 of the ring will flash in yellow twice, then the entire ring will illuminate in yellow.
2	"Hydrating" Mode (Yellow LED)	Quick press the ON/OFF button to turn on the device. After 1 beep, 1/3 of the ring will flash in yellow twice, then the entire ring will illuminate in yellow to show "Hydrating Mode", the device will vibrate and begin to warm.

3	"Brightening" Mode (Red LED)	Whilst in "Hydrating Mode", quickly press the ON/OFF button again to move to "Brightening Mode", you will hear 1 beep, then the device will change to a red light with vibration and warming.
4	"Soothing" Mode (Blue LED)	Whilst in "Brightening Mode", quickly press the ON/OFF button again to move to "Soothing Mode", you will hear 1 beep, then the device will change to a blue light with vibration and warming.
5	OFF	<p>To turn off, press and hold ON/OFF button, you will hear a beep and the product will stop vibrating.</p> <p>The device will automatically shut off after 6 minutes of use in all modes, it will beep and stop vibrating.</p>
6	Low battery	<p>When the battery power is lower than 10%, the white power indicator will flash 3 times, you will be able to use the device once again before charging it.</p> <p>When the battery power is 0%, the white power indicator will flash 3 times and then turn off. You will not be able to use the device again before charging it.</p>

7	Charging	The white power indicator flashing.
8	Fully Charged	The white power indicator will be constantly on.

Specification

No.	Item	Unit	Spec.
1	Nominal Voltage	V	3.7V
2	Operating current	mA	“Hydrating” Mode:40mA-1500mA
			“Brightening” Mode:40mA-1500mA
			“Soothing” Mode:40mA-120mA
3	Charging Time	H	<5H
4	Waterproof	/	IPX7
5	Working Time	Min	6Min

How To Use

The vibrating Miracle Face Mask Device is made from hygienic silicone and innovative photon thermo-therapy technology, gradually warming up to disperse and deeply penetrate the goodness of your favourite facial sheet mask ingredients into the skin, combined with yellow, red and blue LED light which helps to target a variety of skin care concerns, promoting a glowing complexion.


► “Hydrating” Mode:



1. Please cleanse your face before use.



2. Fully open the mask to cover the face, flatten the folds and keep the mask smooth (use a hydrating and moisturising mask).

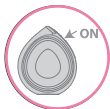

3. Power on and enter "Soothing" mode.



4. Move slowly and massage upwards (Jaw-Cheeks-Nose-Forehead) with the massage plate, avoiding the skin around the eyes.



5. The device will automatically shut down after 6 minutes.



6. Rinse your face with clean, warm water, using a cold faceloth after to close your pores. Follow with your usual moisturiser.



7. Rinse the device with clean water, wipe up and keep it covered.


► “Brightening” Mode:



1. Please cleanse your face before use.



2. Fully open the mask to cover the face, flatten the folds and keep the mask smooth (use a rejuvenating mask).

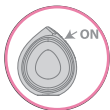

3. Power on and enter "Brightening" mode.



4. Move slowly and massage upwards (Jaw-Cheeks-Nose-Forehead) with the massage plate, avoiding the skin around the eyes.



5. The device will automatically shut down after 6 minutes.



6. Rinse your face with clean, warm water, using a cold faceloth after to close your pores. Follow with your usual moisturiser.



7. Rinse the device with clean water, wipe up and keep it covered.

► “Soothing” Mode:



1. Please cleanse your face before use.



2. Fully open the mask to cover the face, flatten the folds and keep the mask smooth (use a soothing and repairing mask).



3. Power on and enter "Soothing" mode.



4. Move slowly and massage upwards (Jaw-Cheeks-Nose-Forehead) with the massage plate, avoiding the skin around the eyes.



5. The device will automatically shut down after 6 minutes.



6. Rinse your face with clean, warm water, using a cold faceloth after to close your pores.
Follow with your usual moisturiser.



7. Rinse the device with clean water, wipe up and keep it covered.

Tips:

The "▼" mark on the transparent fixed ring indicates the location of the DC charging port.

The device is designed with three modes and it will automatically memorize the last mode used before the last power off.

The waterproof grade of this machine is IPX7.

In the "Hydrating" and "Brightening" modes, the nursing head of the device enters a warm state, and the air inside the product will also generate heat expansion. It is normal for silicone on the back to be slightly bulged.

When the product is not in use, the silicone will go back to its original state.

Warnings

1. Keep out of reach of children.
2. Do not use this device on sunburnt, broken or irritated skin, any swelling or if you have had skin related surgery.
3. Remove watches, glasses, false eyelashes, contact lenses and other solid objects before using.
2. Do not use this appliance if the device is damaged.
4. Do not disassemble or modify this device.
9. Stop use immediately if experiencing any irritation, pain or discomfort.
10. Do not Use this appliance on any infected areas of skin such as rashes or open wounds.
11. Do not use this device during pregnancy.
12. Avoid contact with the eye area.
13. Do not use the device in the same area for too long.
14. Do not share this product with any other person.

Warnings

- ⊘ Do not leave the device unattended while it is working.


ON



- ⚠ Do not use in high temperature environment



- ⊘ Do not disassemble the device by yourself

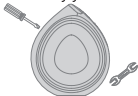

- ⊘ Do not use it if silicone damaged



- ⊘ If the product is heavily hit or dropped to the ground, please check for damage before continuing to use.



Product Care

CLEAN

After each use, rinse the device with clean water immediately, and then dry it up.

STORAGE

Do not immerse the device in water for a long time;

Keep it dry and clean;

Keep it away from high temperature;

Keep it out of reach of children

FAQ's

Q	A	Solution
The device falls to turn on	1.Indicator does not light up (battery or board failure) 2.The white indicator flash (Low battery)	1.Please contact us for help if the device still fails to work after fully charged 2.Restart the device after charging it for 5 hours
The device fails to charging	1.The USB cable is damaged 2.The USB interface is not connected well 3.The battery has reached the end of its life	Please contact us for help if the device still fails to work after replacing the USB cable/power plug.
The device stops working suddenly	1.The device has run out of power 2.There is something wrong with the device	1.Charge for 5 hours and restart to use; 2.Please contact us for help.
NTC failure problem	The mode indicator flashes when the device power on	Temperature sensor failure, return to the factory for repair.

Warm tips: If the product has not been used for a long time (6 months or more), it is recommended to use it after full charge.

SUPPORT

Please feel free to contact us at:



+44 (0) 1483 937967



hello@love-my-skin.co.uk



Love My Skin (Kinetik Medical Devices)

Unit 3, Perrywood Business Park

Honeycrook Lane, Salfords Redhill

RH1 5DZ. UK



LOVE MY SKIN



Symbol for "Manufacture Date".



Symbol for "Environment Protection" – Waste electrical products should not be disposed of with household waste. Please recycle where facilities exist. Check with your local Authority or retailer for recycling advice.



Symbol for "Complies with EU requirements".



Symbol for "Complies with UK requirements".



Follow this instruction manual for your safety.



Model Reference.

Made in China

LMSSM