LOVE MY HAIR

HAIR MASSAGE COMB

to stimulate hair growth

LMS-HMC

Contents

- 3 Introduction
- 4 Product Specifications
- 4 Getting to Know
- 4 Features
- 5 How To Guide
- 6 Keeping it Clean
- 6 Stay Powered
- 7 FAQs
- 8-11 Warnings
- 12 Here to help

Introduction

The Love My Hair, Hair Massage Comb works with 650nm of Red Light Therapy to stimulate the hair follicles by increasing blood circulation and adjust endocrine, used to send the oxygen and nutrition needed to the hair follicle.

This device is also powered with sonic vibration that activates the hair follicle, massages the scalp and relieves fatigue to promote a healthy scalp.


Product Specifications

NAME - Hair Massage Comb

Model Number - LMS-HMC

Power - 5V

Getting to know your Hair Massage Comb

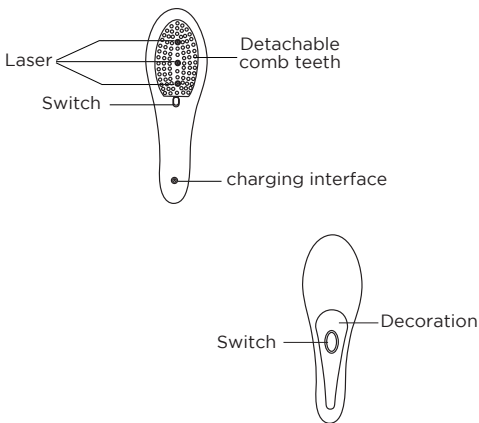Features

- Portable
- USB Charger
- 650nm Red Light Therapy

How To Guide

1. Wash you hair and dry thoroughly
2. Do not apply any product to the hair
3. Press the switch on the back of the appliance lightly, the Red Light and low frequency vibration will begin
4. The vibration intensity can be adjusted by gently pressing the same button again
5. Start from the left ear, moving very slowly, brush the hair from the forehead hairline backwards
6. Repeat across the head to the right ear, and back again
7. Start from the left ear, moving very slowly, brush the hair from one side to the other ear
8. Repeat across the head to the right ear, and back again
9. Turn off the device after use

TOP TIP - The appliance has an auto switch-off at 5 minutes. Use this as a time for stages 5 & 6. Switch the device back on and for a further 5 minutes, use as a timer for stages 7 & 8.

Repeat these movements for 5-10 minutes, 3-4 times per week

Keeping it clean

1. Make sure the appliance is turned off before you begin to clean it
2. Do not immerse the device into water. The detachable parts can be removed for a thorough clean
3. Clean the surface of the device with a dry cloth or with a cloth that is slightly moist
4. Never use lubricants, detergents, polishing agents or air fresheners on this appliance
5. The removable head should be cleaned with a wet cloth
6. Allow to dry before using

Stay powered

1. Connect the small head to the appliance and the USB to a suitable adaptor before plugging into the wall socket
2. When charging has commenced the appliance will have a red light on
3. When charging has completed the appliance will change from a red light to green
4. The charging time is 2 hours
5. Never use the appliance when it is charging

FAQs

Why is my device not working?

Check the appliance is fully charged

Why is my appliance shutting down?

Check the appliance is fully charged. If the appliance has been on for 5 minutes, it will automatically switch-off. Simply press the 'ON' button again and continue use for another 5 minutes.

Can I use this if I have sensitive skin or a skin condition?

Whilst the product is very gentle, if you have a pre-existing condition, we do recommend you consult your doctor or dermatologist before using.

Warnings

- KEEP OUT OF REACH OF CHILDREN
- Do not look directly into the Red Light
- Do not attempt to disassemble or repair the device
- Do not place near fire or in an area of high temperature
- Ensure the appliance is cleaned thoroughly after each use
- If not using for long periods of time, it is recommended to still fully charge the appliance every 6 months at a minimum
- Do not submerge the appliance directly into water
- If appliance begins to feel uncomfortable, immediately stop use

Not suitable for people as follows

- Those suffering with skin allergies, skin diseases, eczema or other skin conditions
- Do not use if pregnant or breast-feeding
- Do not use if you have recently suffered sunburn or high exposure to the sun
- Do not use if you have a pacemaker or implants without first consulting your doctor
- Do not use if you are sensitive to light or use photosensitizing drugs without first consulting your doctor
- Do not use if you have diabetes or high blood pressure without fist consulting your doctor
- If you are on medication, always consult your doctor first

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 hello@love-my-skin.co.uk

 www.love-my-hair.co.uk

 Love My Hair (Kinetik Medical Devices Ltd)
Unit 3, Perrywood Business Park
Honeycrock Lane
Salfords, Redhill
RH1 5DZ
UK

LOVE MY HAIR

 Symbol for “Manufacture Date”

 Symbol for “Environment Protection” - Waste electrical products should not be disposed of with household waste. Please recycle where facilities exist. Check with your local authority or retailer for recycling advice.

 Symbol for “Complies with EU Requirements”.

 Symbol for “Complies with UK Requirements”.

 Follow this instruction manual for your safety.

 Symbol for “No use in bathroom”.

 Model Reference.

Made in China
LMS-HMC